



IRRIGATION SET
ASPIRATION SET



VECTEC
SINGLE USE
LAPAROSCOPIC
DEVICES



VECTEC




REF : RAL 3000

GRAVITY BASED IRRIGATION/ASPIRATION SET to fit general and laparoscopic surgery

PRODUCTS

Reference	Designation	Diameter	Length	Packaging
CAN250 5	irrigation/aspiration set	5 mm	330 mm	10 units per box
CAN250 5P	irrigation/aspiration set with anti-drip tip	5 mm	330 mm	10 units per box
CAN250 5W	irrigation/aspiration set with adaptor	5 mm	330 mm	10 units per box
CAN250 5GL	irrigation/aspiration set	5 mm	450 mm	10 units per box
CAN250 5DP	irrigation/aspiration set with double spike	5 mm	330 mm	10 units per box
CAN250 10	irrigation/aspiration set	10 mm	330 mm	10 units per box
CAN 5LS	Suction irrigation aspiration with cannula	5 mm	330 mm	10 units per box

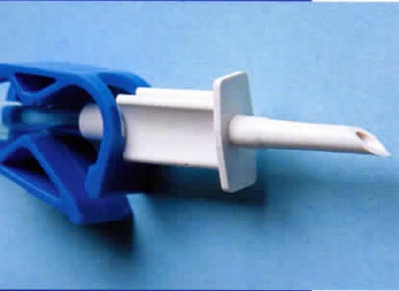

Spike connector



Adaptor



Ø5 mm/10 mm cannula

ACCESSORIES

Reference	Designation	Diameter	Length	Packaging
CAN 10	Cannula only	10 mm	330 mm	10 units per box
RAL 3000	Extension suction tubing	7 mm	3 m	25 units per box

The single use suction irrigation set consists in :

- An atraumatic interchangeable 5 mm or 10 mm sheath with distal lateral holes.
- An ergonomic handle fitted with double trumpet-valves with colour coded plunger, also suitable for biopsy
- A transparent and flexible 2,5 meter long double tubing set
- A single spike connection for irrigation bags and an on-and-off clamp
- A soft universal connector fitting to any suction canister

MD Class IIa
Sterilized by ETO to
NF EN ISO 11135-1 standards
See instructions for use
Manufactured by Promepla

CE
0459

G
MED NF EN ISO 13485



REF : CAN 10

Distributed by :

VECTEC

Bioparc - Vichy - Hauterive
03270 Hauterive
FRANCE

Phone : +33 (0)4 70 58 41 70

Fax : +33 (0)4 70 58 41 79

Email : vectec@vectec.fr

http://www.vectec.fr